

# LanthaScreen™ TR-FRET Assay Performance with Thermo Scientific Varioskan Flash Spectral Scanning Multimode Reader

Jorma Lampinen, Thermo Fisher Scientific, Vantaa, Finland and Kevin Lowitz, Invitrogen Corporation, Carlsbad, California, USA



LanthaScreen™ Certified



invitrogen

## Summary

LanthaScreen TR-FRET assay technology from Invitrogen is a widely used application in biotech and pharmaceutical research. This paper describes the performance of the Varioskan® Flash multimode reader using LanthaScreen assays. The results show that Varioskan Flash is compatible with LanthaScreen assays, and Varioskan Flash is an excellent spectral scanning multimode reader for TR-FRET based applications.

## Introduction

LanthaScreen technology uses time-resolved fluorescence resonance energy transfer (TR-FRET) to study multiple target classes including protein kinases, nuclear hormone receptors, as well as proteases and ubiquitinated proteins. LanthaScreen TR-FRET works on the principles that when suitable pairs of fluorophores are in close proximity of one another,

excitation of the terbium chelate donor fluorophore results in energy transfer to the fluorescein or GFP acceptor fluorophore (Figure 1). The benefit of this assay is the long fluorescence

lifetime of the terbium lanthanide donor fluorochrome excited at 340 nm. This allows ratiometric measurement of acceptor fluorescence emission (FRET signal) at 520 nm over the 495 nm donor emission long after the background fluorescence has dissipated (Figure 2). The time-resolved measurement reduces assay interference (e.g. fluorescent compounds) and increases data quality.

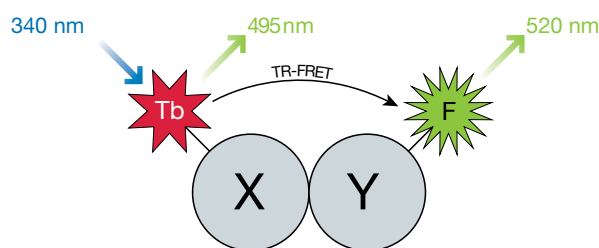


Figure 1. Principle of LanthaScreen TR-FRET detection

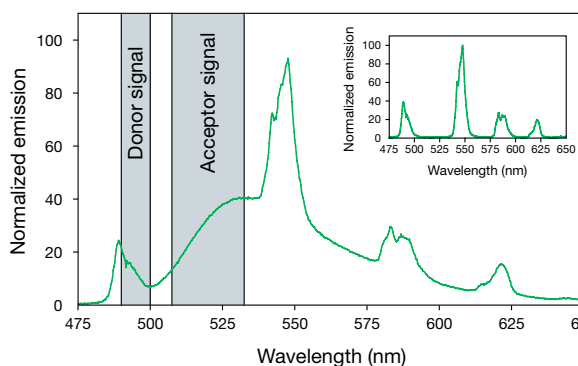


Figure 2. Fluorescence Spectra of Terbium and Fluorescein (FRET) emission. The time-resolved spectra above illustrate energy transfer occurring when terbium (donor) and fluorescein (acceptor) are brought into proximity via biomolecular interactions. The inset shows the time-resolved spectra in the absence of energy transfer.

## Materials and methods

The performance of Varioskan Flash was tested using a LanthaScreen Tb-anti-GST Antibody Kit (catalogue number PV4216) obtained from Invitrogen (Carlsbad, CA). Solid white 384-well microplates were purchased from Corning (Corning, NY). The assay was performed according to the kit instructions.

Briefly, the fluorescein labeled GST-MBP positive control provided was serially diluted from 200nM stock solution and the Tb-anti-GST antibody was diluted to 4 nM concentration using the LanthaScreen TR-FRET dilution buffer provided in the kit. For reactions in standard white 384-well plates 20,  $\mu\text{L}$  of Tb-anti-GST antibody (2nM final concentration) and 20  $\mu\text{L}$  of the fluorescein-GST-MBP positive control dilution were aliquoted in the plate with four replicates. When low volume white 384-well plates were used, aliquots of 10  $\mu\text{L}$  of Tb-anti-GST antibody and 10  $\mu\text{L}$  of the fluorescein-GST-MBP positive control dilution were added into the plate using the same number of replicates. The reactions were allowed to incubate for 60 minutes at room temperature and then measured with the Varioskan Flash controlled with Thermo Scientific SkanIt Software (Thermo Fisher Scientific, Finland). The measurement parameters programmed into the SkanIt<sup>®</sup> Software to perform time-resolved fluorescence dual emission measurement are shown in Table I. Required dual emission measurement was performed using the software's multi-wavelength measurement option to measure both donor and acceptor emissions easily from the sample without any additional time differences. The structure of the measurement protocol and multi-wavelength measurement step used is shown in Figure 3.



Figure 3. SkanIt Software measurement protocol and parameters for LanthaScreen TR-FRET assays.

	Measurement 1 (Tb-donor)	Measurement 2 (Fluorescein acceptor)
Excitation wavelength (nm)	332	332
Emission wavelength (nm)	488	518
Excitation bandwidth (nm)	12	12
Emission bandwidth (nm)	12	12
TRF delay time (us)	100	100
TRF integration time (us)	200	200
Measurement time (ms)	1000 (100 flashes)	1000 (100 flashes)

Table I. Varioskan Flash measurement parameters for LanthaScreen assays

After the measurement, the data from both standard and low volume 384-well plate assays was transferred to Microsoft Excel<sup>®</sup> for further calculations. The TR-FRET ratio between acceptor and donor emission signals were calculated for each positive control sample and calibration curves were drawn. Then, theoretical assay sensitivities for both plate formats using standard IUPAC  $3 \times \text{SD}$  principle and Z-prime values for each positive control concentration were calculated.

## Results

Calibration curves of the assays with 384- and low volume 384-well plates are shown in Figure 4, while Figure 5, shows Z-prime values against positive control concentrations.

According to these results, both standard and low volume white 384-well plates perform well in LanthaScreen assays with the Varioskan Flash. Measured calibration curves show almost

identical behavior with both plate types and very similar theoretical assay sensitivities are obtained with both plate formats.

When assay robustness is estimated based on Z-prime values, similar behavior is also noticed between standard and low volume 384-well plates. Z-prime values are well over required limit of 0.5 with all sample concentrations over 1 nM.

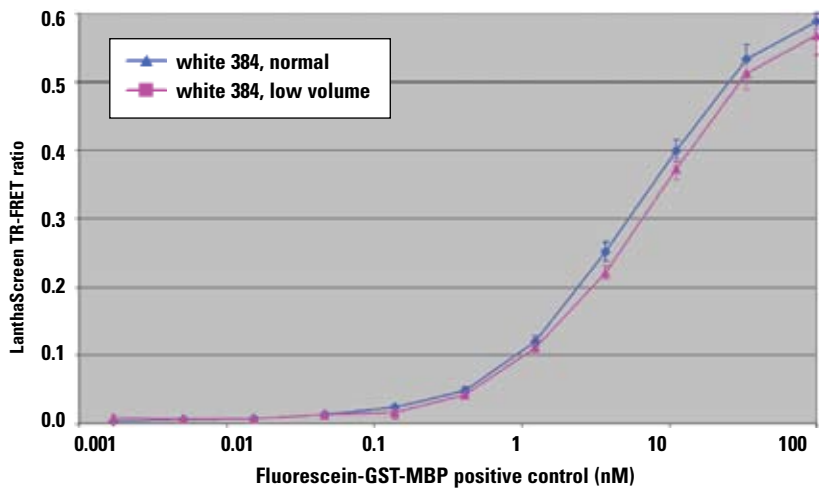


Figure 4. Calibration curves of LanthaScreen assays with standard and low volume white 384-well plates

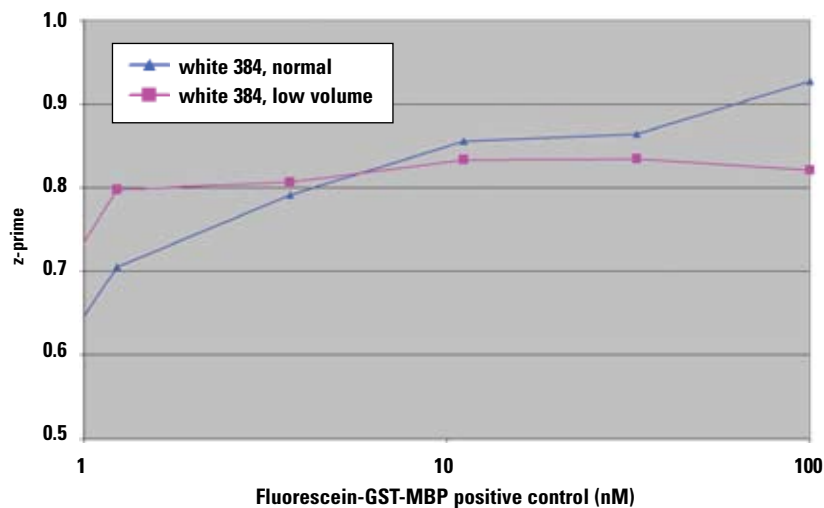


Figure 5. Z-prime values of LanthaScreen assays with Varioskan Flash

## Conclusions

These results show that the Varioskan Flash fulfills all requirements for LanthaScreen assays and can be used successfully to measure TR-FRET based assays with good sensitivity. Furthermore, measuring TR-FRET assays that require dual emission measurement is easy to program and can be done with only one measurement step using the software's multi-wavelength feature. LanthaScreen assays can also be performed with Varioskan Flash using low total assay volume and low volume 384-well plates that makes this Varioskan Flash/LanthaScreen combination especially useful for screening purposes where low sample consumption is exceptionally important.

## Further information

For further information about the Varioskan Flash spectral scanning multimode reader or Invitrogen LanthaScreen technology, please refer to the following web pages:

[www.thermo.com/readingroom](http://www.thermo.com/readingroom)

[www.thermo.com/varioskan](http://www.thermo.com/varioskan)

[www.invitrogen.com/lanthascreen](http://www.invitrogen.com/lanthascreen)

*In addition to these offices, Thermo Fisher Scientific maintains a network of representative organizations throughout the world.*

### North America:

**USA / Canada**  
+1 800 522 7763

### Europe:

**Austria**  
+43 1 801 40 0

**Belgium**  
+32 2 482 30 30

**France**  
+33 2 2803 2180

**Germany national toll free**  
08001-536 376

**Germany international**  
+49 6184 90 6940

**Italy**  
+39 02 02 95059 448

**Netherlands**  
+31 76 571 4440

**Nordic countries**  
+358 9 329 100

**Russia/CIS**  
+7 (495) 739 76 41

**Spain/Portugal**  
+34 93 223 09 18

**Switzerland**  
+41 44 454 12 12

**UK/Ireland**  
+44 870 609 9203

### Asia:

**China**  
+86 21 6865 4588 or  
+86 10 8419 3588

**India toll free**  
1800 22 8374

**India**  
+91 22 6716 2200

**Japan**  
+81 45 453 9220

**Other Asian countries**  
+852 2885 4613

**Countries not listed:**  
+49 6184 90 6940 or  
+33 2 2803 2180

[www.thermo.com/readingroom](http://www.thermo.com/readingroom)  
[www.thermo.com/mpi](http://www.thermo.com/mpi)

© 2008 Thermo Fisher Scientific Inc. All rights reserved. LanthaScreen is a trademark of Invitrogen Corporation. Excel is a registered trademark of Microsoft Corp. All other trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.