


Shandon EZ Cytofunnel[®]


Bench Stand – Kit A78710022

Carefully position the EZ Cytofunnel[®], with the slide fitted, into a slot.

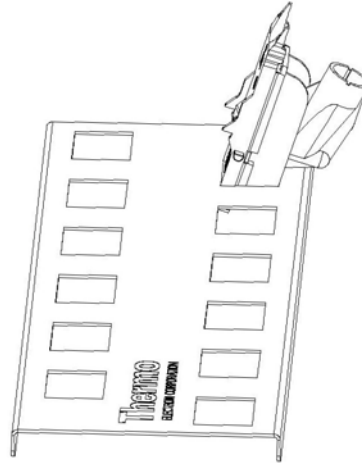

Pipette the appropriate volume of fluid into the chamber.

Carefully remove the EZ Cytofunnel[®] and insert into the Cytospin[®] Head.

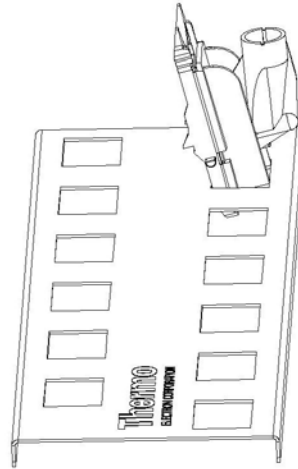

When handling the EZ Cytofunnel[®] with fluid in the chamber, do not allow the fluid to contact the slide.

Thermo

ELECTRON CORPORATION

Thermo Electron Corporation, Anatomical Pathology

93-96 Chadwick Road, Astmoor, Runcorn, Cheshire, WA7 1PR, United Kingdom

Tel: +44 (0) 1928 566611; Fax: +44(0) 1928 565845; www.thermo.com/shandon